



# *EMS Vacuum Immobilisation Catalogue 2021*



# **RedVac EMS Vacuum Immobilisation**

## **Content Page**

**RedVac made in Austria 3**

**The RedVac System 4-5**

**Mattresses 6-10**

**Head Immobilisation 11**

**Extremity Splints 12-16**

**Extrication Device 17**

**Pediatric 18-19**

**Accessories 20**

**Pumps & Bags 21**

**Thermo Rescue Sack 22**

**Helicopter Rescue Sack 23**



# Developed and manufactured in Austria

Kohlbrat & Bunz products are developed and manufactured 900 meters above sea level in the middle of the Austrian Alps in its own facilities by its experienced and meticulously working staff.

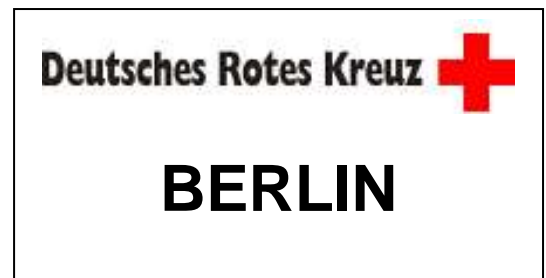
People living in the Alps know very well that products used under the extreme conditions of alpine rescue have to work perfectly.

To meet these rigorous demands our products are developed in close cooperation with the mountain rescue teams and other demanding emergency services.

This cooperation has brought about a selection of the best materials and optimal product design leading to products of outstanding performance and quality. To insure the best possible quality Kohlbrat & Bunz work in accordance with EN ISO 9001 & EN ISO 13485.

## Worldwide in use

We are very proud of the confidence shown to us and our products by many well-known and demanding emergency services worldwide and are happy to show some examples of their beautiful logos that they have asked us to put on the RedVac products we have made for them.



# The RedVac System

for  
**PERFECT HYGIENE & EASE OF USE**



- **The RedVac CLICK-System**  
*for attachment of straps and handles*

*The unique CLICK-System makes removal of straps and handles for repositioning, replacement, cleaning and disinfection very easy!*



- **The RedVac Multi-Chamber-Systems**  
*- provides unmatched EASE OF USE*

*The Multi-Chamber-Systems keeps the internal beads evenly distributed eliminating the need for difficult and time consuming manual spreading of the beads. This enables fast, safe and easy application in all situations.*

**The RedVac Multi-Chamber-System!**

*The internal sack system keeps the beads evenly distributed in up to 120 chambers.*



**The RedVac DOT-Chamber-System!**

*The layers of the outer cover are welded together dot by dot creating chambers between the dots keeping the internal beads in their correct position, evenly distributed.*



- **Optional attachment systems for straps and handles**

**Cleats**



**Hook & Loop**



**Red Stripes**



- **Robust & durable exterior Material**

To achieve outstanding durability, abrasion resistance, tightness, easy cleaning and disinfection, comfort and safety for the patient special materials were chosen. It makes it possible to use the RedVac products in a wide temperature range, from -40°C up to +80°C.



- **Restraint Straps, made of washable textile or wipe-clean coated fabric**

Colour coded patient restraint straps, adjustable from both sides, easy to clean, reposition or replace.



- **Vacuum Valve, screw-on replaceable**

The valve is robust and easy to use, even when wearing gloves. It closes automatically when suction stops preventing loss of vacuum if slipping with the pump. Comes with suitable adapter for most vacuum pumps.



- **Vacuum Valve, welded-on non-replaceable**

The valve is robust and easy to use, even when wearing gloves. It closes automatically when suction stops preventing loss of vacuum if slipping with the pump.



- **Label System**

An indestructible label with individual serial number and date of manufacture is welded onto every product. Optional labels with property indicator and distributor logo are available.



- **X-ray, CT and MRI compatible**



- **CE**

All RedVac products are CE marked and conform with the European directive 93/42/EEC for medical products class 1.




- **ISO9001 - ISO13485 certified by**




# RedVac EMS Vacuum-Mattresses




## Body shape

	<b>VM01X02</b>	Cleats with webbing restraint straps
	<b>VM01CX02</b>	CLICK with webbing restraint straps
	<b>VM01CHX02</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 40-88-50cm Length: 200cm Weight: 7kg Carry handles: 8 Patient restraint straps: 4	


## Body shape

	<b>VM03X03</b>	Cleats with webbing restraint straps
	<b>VM03CX03</b>	CLICK with webbing restraint straps
	<b>VM03CHX03</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 50-88-50cm Length: 200cm Weight: 7kg Carry handles: 8 Patient restraint straps: 4	


## Body shape. Slits for gurney safety belts in neck area.

	<b>VM03X05</b>	Cleats with webbing restraint straps
	<b>VM03CX05</b>	CLICK with webbing restraint straps
	<b>VM03CHX05</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 50-98-50cm Length: 200cm Weight: 7kg Carry handles: 8 Patient restraint straps: 4	

## Body shape. Slits for gurney safety belts in neck area. Glass fiber rod for cervical protection.


	<b>VM03X07</b>	Cleats with webbing restraint straps
	<b>VM03CX07</b>	CLICK with webbing restraint straps
	<b>VM03CHX07</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 50-88-50cm Length: 200cm Weight: 7kg Carry handles: 8 Patient restraint straps: 4	

## Body shape - compact size for helicopter transportation


	<b>VM08X06</b>	Cleats with webbing restraint straps
	<b>VM08CX06</b>	CLICK with webbing restraint straps
	<b>VM08CHX06</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 40-78-40cm Length: 195cm Weight: 6kg Carry handles: 8 Patient restraint straps: 4	
	<b>VM09X01</b>	Cleats with webbing restraint straps
	<b>VM09CX01</b>	CLICK with webbing restraint straps
	<b>VM09CHX01</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 40-60-40cm Length: 195cm Weight: 5kg Carry handles: 8 Patient restraint straps: 4	

## Body shape with head fixation.

### Slits for gurney safety belts in neck area. Glass fiber rod for cervical protection.

	<b>VM200X06</b>	Cleats with webbing restraint straps
	<b>VM200CX06</b>	CLICK with webbing restraint straps
	<b>VM200CHX06</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 40-88-60cm Length: 200cm Weight: 7kg Carry handles: 8 Patient restraint straps: 6	
	<b>VM201X10</b>	Cleats with webbing restraint straps
	<b>VM201CX10</b>	CLICK with webbing restraint straps
	<b>VM201CHX10</b>	CLICK with wipe-clean restraint straps
	Width head/mid/foot: 40-100-60cm Length: 210cm Weight: 8kg Carry handles: 8 Patient restraint straps: 6	

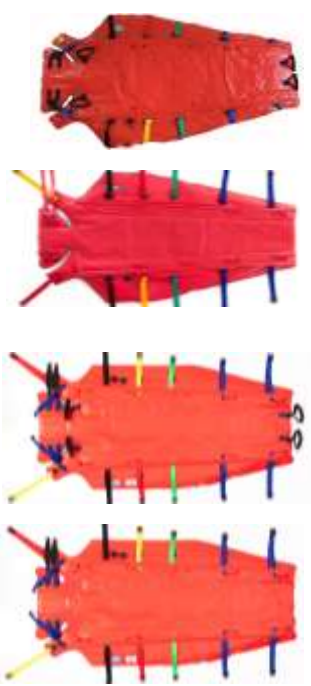
**Head fixation**

	<p><b>VM500X03</b> Cleats with webbing restraint straps  <b>VM500CX03</b> CLICK with webbing restraint straps  <b>VM500CHX03</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 60-88-88cm Length: 200cm Weight: 7,5kg          Carry handles: 8 Patient restraint straps: 5</p>
	<p><b>VM501X03</b> Cleats with webbing restraint straps  <b>VM501CX03</b> CLICK with webbing restraint straps  <b>VM501CHX03</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 60-100-100cm Length: 210cm Weight: 8,5kg          Carry handles: 8 Patient restraint straps: 5</p>
	<p><b>VM502X01</b> Cleats with webbing restraint straps  <b>VM502CX01</b> CLICK with webbing restraint straps  <b>VM502CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 60-80-80cm Length: 200cm Weight: 7kg          Carry handles: 8 Patient restraint straps: 5</p>



**Head fixation and slits in neck area for gurney safety belts. Glassfiber rod for cervical protection.**

	<p><b>VM5000X01</b> Cleats with webbing restraint straps  <b>VM5000CX01</b> CLICK with webbing restraint straps  <b>VM5000CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 65-88-88cm Length: 200cm Weight: 7kg          Carry handles: 8 Patient restraint straps: 5</p>
	<p><b>VM5001X01</b> Cleats with webbing restraint straps  <b>VM5001CX01</b> CLICK with webbing restraint straps  <b>VM5001CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 65-100-100cm Length: 210cm Weight: 8kg          Carry handles: 12 Patient restraint straps: 5</p>
	<p><b>VM5002X01</b> Cleats with webbing restraint straps  <b>VM5002CX01</b> CLICK with webbing restraint straps  <b>VM5002CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 65-80-80cm Length: 200cm Weight: 6kg          Carry handles: 8 Patient restraint straps: 5</p>


**Body shape with head fixation. Slits for gurney safety belts. Glass fiber rod for cervical protection. Additional outer protective layer with 10 integrated handles.**

	<p><b>VM7000X03</b> Cleats with webbing restraint straps  <b>VM7000CX03</b> CLICK with webbing restraint straps  <b>VM7000CHX03</b> CLICK with wipe-clean restraint straps</p> <p>Width: 130cm Length: 205cm Weight: 10kg          Carry handles non-removable: 10 Carry handles removable: 4          Patient restraint straps: 5</p>
	<p><b>VM7000X04</b> Cleats with webbing restraint straps  <b>VM7000CX04</b> CLICK with webbing restraint straps  <b>VM7000CHX04</b> CLICK with wipe-clean restraint straps</p> <p>Width: 130cm Length: 205cm Weight: 10kg          Carry handles non-removable: 10 Patient restraint straps: 5</p>
	<p><b>VM7001X01</b> Cleats with webbing restraint straps  <b>VM7001CX01</b> CLICK with webbing restraint straps  <b>VM7001CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width: 108cm Length: 205cm Weight: 9,5kg          Carry handles non-removable: 10 Carry handles removable: 4          Patient restraint straps: 5</p>
	<p><b>VM7001X04</b> Cleats with webbing restraint straps  <b>VM7001CX04</b> CLICK with webbing restraint straps  <b>VM7001CHX04</b> CLICK with wipe-clean restraint straps</p> <p>Width: 110cm Length: 205cm Weight: 9,5kg          Carry handles non-removable: 10 Patient restraint straps: 5</p>

**Head fixation. Slits for gurney safety belts in neck area. Glassfiber rod for cervical protection**

	<p><b>VM11X04</b> Cleats with webbing restraint straps  <b>VM11CX04</b> CLICK with webbing restraint straps  <b>VM11CHX04</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 80-88-80cm Length: 200cm Weight: 7kg          Carry handles: 8 Restraint straps: 4</p>
	<p><b>VM1101X01</b> Cleats with webbing restraint straps  <b>VM1101CX01</b> CLICK with webbing restraint straps  <b>VM1101CHX04</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 80-100-80cm Length: 210cm Weight: 8kg          Carry handles: 12 Restraint straps: 4</p>
	<p><b>VM1111X01</b> Cleats with webbing restraint straps  <b>VM1111CX01</b> CLICK with webbing restraint straps  <b>VM1111CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 80-108-80cm Length: 210cm Weight: 8,5kg          Carry handles: 12 Restraint straps: 4</p>
	<p><b>VM1112X01</b> Cleats with webbing restraint straps  <b>VM1112CX01</b> CLICK with webbing restraint straps  <b>VM1112CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width head/mid/foot: 80-108-80cm Length: 220cm Weight: 8,5kg          Carry handles: 8 Restraint straps: 4</p>

**Rectangular shape. Slits for gurney safety belts in the neck area Glassfiber rod for cervical protection.**


	<p><b>VM11620X01</b> Cleats with webbing restraint straps  <b>VM11620CX01</b> CLICK with webbing restraint straps  <b>VM11620CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width: 60cm Length: 200cm Weight: 5kg          Carry handles: 8 Restraint straps: 4</p>
	<p><b>VM11720X01</b> Cleats with webbing restraint straps  <b>VM11720CX01</b> CLICK with webbing restraint straps  <b>VM11720CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width: 70cm Length: 200cm Weight: 6kg          Carry handles: 8 Restraint straps: 4</p>
	<p><b>VM11820X01</b> Cleats with webbing restraint straps  <b>VM11820CX01</b> CLICK with webbing restraint straps  <b>VM11820CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width: 80cm Length: 200cm Weight: 7kg          Carry handles: 8 Restraint straps: 4</p>
	<p><b>VM11920X01</b> Cleats with webbing restraint straps  <b>VM11920CX01</b> CLICK with webbing restraint straps  <b>VM11920CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width: 90cm Length: 200cm Weight: 7,5kg          Carry handles: 8 Restraint straps: 4</p>
	<p><b>VM1201X01</b> Cleats with webbing restraint straps  <b>VM1201CX01</b> CLICK with webbing restraint straps  <b>VM1201CHX01</b> CLICK with wipe-clean restraint straps</p> <p>Width: 100cm Length: 210cm Weight: 8,5kg          Carry handles: 8 Restraint straps: 4</p>



### Tapering shape at head and foot area

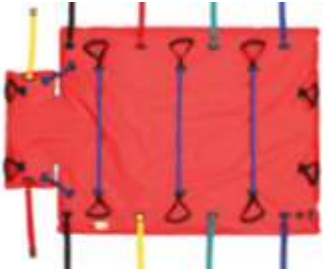

	<b>VM10X02</b> Cleats with webbing restraint straps <b>VM10CX02</b> CLICK with webbing restraint straps <b>VM10CHX02</b> CLICK with wipe-clean restraint straps Width head/mid/foot: 80-88-80cm Length: 200cm Weight: 7kg Carry handles: 8 Restraint straps: 4
	<b>VM8021X01</b> Cleats with webbing restraint straps <b>VM8021CX01</b> CLICK with webbing restraint straps <b>VM8021CHX01</b> CLICK with wipe-clean restraint straps Width head/mid/foot: 80-100-80cm Length: 210cm Weight: 8,5kg Carry handles: 12 Restraint straps: 4
	<b>VM8121X01</b> Cleats with webbing restraint straps <b>VM8121CX01</b> CLICK with webbing restraint straps <b>VM8121CHX01</b> CLICK with wipe-clean restraint straps Width head/mid/foot: 80-108-80cm Length: 210cm Weight: 9kg Carry handles: 12 Restraint straps: 4

### Rectangular shape

	<b>VM620X01</b> Cleats with webbing restraint straps <b>VM620CX01</b> CLICK with webbing restraint straps <b>VM620CHX01</b> CLICK with wipe-clean restraint straps Width: 60cm Length: 200cm Weight: 6kg Carry handles: 8 Restraint straps: 4
	<b>VM720X01</b> Cleats with webbing restraint straps <b>VM720CX01</b> CLICK with webbing restraint straps <b>VM720CHX01</b> CLICK with wipe-clean restraint straps Width: 70cm Length: 200cm Weight: 6,5kg Carry handles: 8 Restraint straps: 4
	<b>VM820X01</b> Cleats with webbing restraint straps <b>VM820CX01</b> CLICK with webbing restraint straps <b>VM820CHX01</b> CLICK with wipe-clean restraint straps Width: 80cm Length: 200cm Weight: 7kg Carry handles: 8 Restraint straps: 4
	<b>VM1021X01</b> Cleats with webbing restraint straps <b>VM1021CX01</b> CLICK with webbing restraint straps <b>VM1021CHX01</b> CLICK with wipe-clean restraint straps Width: 100cm Length: 210cm Weight: 8,5kg Carry handles: 12 Restraint straps: 4
	<b>VM1121X01</b> Cleats with webbing restraint straps <b>VM1121CX01</b> CLICK with webbing restraint straps <b>VM1121CHX01</b> CLICK with wipe-clean restraint straps Width: 108cm Length: 210cm Weight: 9kg Handles: 12 Carry handles: 12 Restraint straps: 4

**For obese patients**

**Slits for gurney safety belts in the neck area Glassfiber rod for cervical protection.**


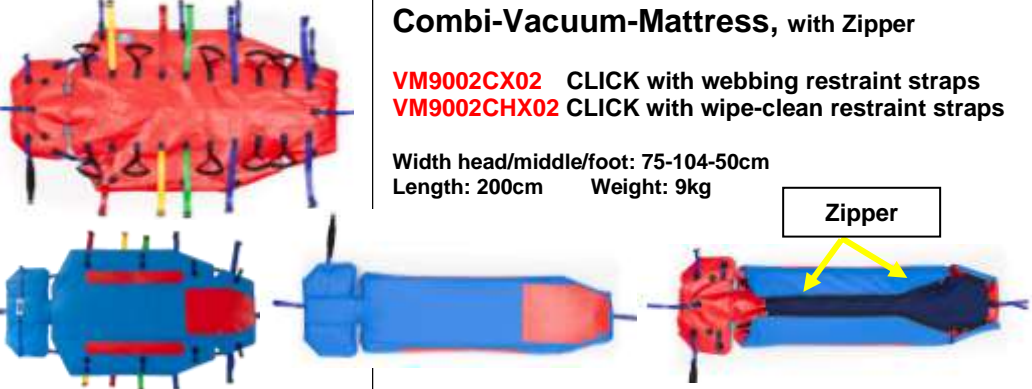

	<p><b>VM1421X01</b>      Cleats with webbing restraint straps</p>
	<p><b>VM1421CX01</b>      CLICK with webbing restraint straps  <b>VM1421CHX01</b>      CLICK with wipe-clean restraint straps</p>

**The RedVac Combi-Vacuum-Mattress**

**2 in 1 - Combines the standard mattress and the vacuum mattress**

**Body shape with head fixation. Slits for gurney safety belts.**

**Glass fiber rod for cervical protection.**

	<p><b>Combi-Vacuum-Mattress, without Zipper</b></p> <p><b>VM9001CX01</b>      CLICK with webbing restraint straps</p> <p><b>VM9001CHX01</b>      CLICK with wipe-clean restraint straps</p> <p>Width head/middle/foot: 75-104-50cm          Length: 200cm      Weight: 9kg</p>
	<p><b>Combi-Vacuum-Mattress, with Zipper</b></p> <p><b>VM9002CX02</b>      CLICK with webbing restraint straps  <b>VM9002CHX02</b>      CLICK with wipe-clean restraint straps</p> <p>Width head/middle/foot: 75-104-50cm          Length: 200cm      Weight: 9kg</p>
	<p><b>Vacuum Pillow</b>  <b>VK200</b></p> <p>Width: 57cm      Length: 39cm      Weight: 0,5kg</p>



## The RedVac Vacuum-Head-Immobiliser - for the Spine Board & the Scoop Stretcher

- Quick and easy positioning of the immobiliser
- No manipulation of head and cervical spine
- Adjusts automatically to any head shape and size
- Ears remain free for unrestricted hearing and medical examination
- Solid and comfortable immobilisation



### Specification

	<p><b>Head Immobiliser</b> <b>VK98X07</b> Width: 38cm Length: 31cm Weight: 0,7kg</p> <p><b>Base plates</b> <b>GP08</b> for scoop stretcher <b>GP11</b> for back board</p>
	<p><b>Head Immobiliser with buckles for attachment</b> <b>VK94X17</b></p> <p>Width: 38cm      Length: 31cm Weight: 0,7kg</p>



## RedVac Vacuum Extrication Device

- the answer to your extrication challenges!
- fits every patient, pregnant – overweight ... no problem!
- comfortable – no respiratory impairment!

### Specification

<p><b>VED24X02</b></p> <p>Width, head: 40cm Width, middle: 106cm Length: 94cm Weight: 2,5kg</p>
-------------------------------------------------------------------------------------------------------------



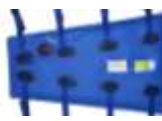
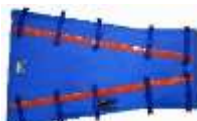
# RedVac EMS Vacuum Extremity Splints



## Optional restraint systems:


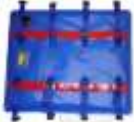






CLICK-System Hook & Loop Red Stripes Cleats







Patient straps:  
Webbing of PP  
Wipe-clean of PVC  
Choice of colours




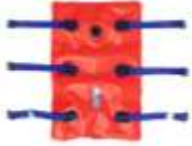

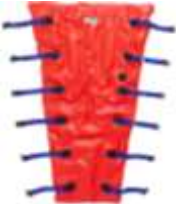


### Angulate design







	<b>Forearm Splints</b> <b>VS60X01</b> Red Stripes with Webbing <b>VS61X01</b> Red Stripes with Webbing Restraint straps: 3 Width: 35cm Length: 50cm Weight: 0,33kg	1-Chamber Multi-Chamber-Sack
	<b>Arm Splints</b> <b>VS62X01</b> Red Stripes with Webbing <b>VS63X03</b> Red Stripes with Webbing Restraint straps: 5 Width top/foot: 43-70cm Length: 56cm Weight: 0,9kg	1-Chamber Multi-Chamber-Sack
	<b>Leg Splints</b> <b>VS64X21</b> Red Stripes with Webbing <b>VS65X21</b> Red Stripes with Webbing Restraint straps: 5 Width top/foot: 64-70cm Length: 103cm Weight: 1,9kg	1-Chamber Multi-Chamber-Sack
	<b>Forearm Splints</b> <b>VS60CX05</b> CLICK with Webbing <b>VS60CHX05</b> CLICK with Wipe-Clean Straps <b>VS61CX05</b> CLICK with Webbing <b>VS61CHX05</b> CLICK with Wipe-Clean Straps <b>VS61CPX05</b> CLICK with Webbing <b>VS61CPHX05</b> CLICK with Wipe-Clean Straps Restraint straps: 3 Width: 35cm Length: 50cm Weight: 0,33kg	1-Chamber 1-Chamber Multi-Chamber-Sack Multi-Chamber-Sack DOT-Chamber DOT-Chamber
	<b>Arm Splints</b> <b>VS62CX05</b> CLICK with Webbing <b>VS62CHX05</b> CLICK with Wipe-Clean Straps <b>VS63CX05</b> CLICK with Webbing <b>VS63CHX05</b> CLICK with Wipe-Clean Straps <b>VS63CPX05</b> CLICK with Webbing <b>VS63CPHX05</b> CLICK with Wipe-Clean Straps Restraint straps: 5 Width top/foot: 43-70cm Length: 56cm Weight: 0,9kg	1-Chamber 1-Chamber Multi-Chamber-Sack Multi-Chamber-Sack DOT-Chamber DOT-Chamber
	<b>Leg Splints</b> <b>VS64CX05</b> CLICK with Webbing <b>VS64CHX05</b> CLICK with Wipe-Clean Straps <b>VS65CX23</b> CLICK with Webbing <b>VS65CHX23</b> CLICK with Wipe-Clean Straps <b>VS65CPX23</b> CLICK with Webbing <b>VS65CPHX23</b> CLICK with Wipe-Clean Straps Restraint straps: 5 Width top/foot: 64-70cm Length: 103cm Weight: 1,9kg	1-Chamber 1-Chamber Multi-Chamber-Sack Multi-Chamber-Sack DOT-Chamber DOT-Chamber
	<b>Forearm Splints</b> <b>VS60X20V</b> Hook & Loop straps <b>VS61X20V</b> Hook & Loop straps <b>VS61PX20V</b> Hook & Loop straps Restraint straps: 3 Width: 35cm Length: 50cm Weight: 0,33kg	1-Chamber Multi-Chamber-Sack DOT-Chamber
	<b>Arm Splints</b> <b>VS62X20V</b> Hook & Loop straps <b>VS63X20V</b> Hook & Loop straps <b>VS63PX20V</b> Hook & Loop straps Restraint straps: 4 Width top/foot: 43-70cm Length: 56cm Weight: 0,9kg Welded-on non-removable valve	1-Chamber Multi-Chamber-Sack DOT-Chamber
	<b>Leg Splints</b> <b>VS64X20V</b> Hook & Loop straps <b>VS65X20V</b> Hook & Loop straps <b>VS65PX20V</b> Hook & Loop straps Restraint straps: 5 Width top/foot: 64-70cm Length: 103cm Weight: 1,9kg Welded-on non-removable valve	1-Chamber Multi-Chamber-Sack DOT-Chamber

	<p><b>Arm splint</b></p> <p><b>VA14X02</b> Red Stripes with Webbing 1-Chamber  <b>VA14X01</b> Red Stripes with Webbing Multi-Chamber-Sack</p> <p>Restraint straps: 4 Width: top/hand: 50-30cm Length: 67cm Weight: 0,9kg</p>
	<p><b>Ankle splint</b></p> <p><b>VK15X02</b> Red Stripes with Webbing 1-Chamber  <b>VK15X01</b> Red Stripes with Webbing Multi-Chamber-Sack</p> <p>Restraint straps: 4 Width: 47cm Length: 57cm Weight: 0,9kg</p>
	<p><b>Leg splint</b></p> <p><b>VL21X02</b> Red Stripes with Webbing 1-Chamber  <b>VL21X01</b> Red Stripes with Webbing Multi-Chamber-Sack</p> <p>Restraint straps: 5 Width: top/hand: 70-40cm Length: 110cm Weight: 1,9kg</p>
	<p><b>Large leg splint</b></p> <p><b>VL18X02</b> Red Stripes with Webbing 1-Chamber  <b>VL18X01</b> Red Stripes with Webbing Multi-Chamber-Sack</p> <p>Restraint straps: 5 Width: top/hand: 70-45cm Length: 135cm Weight: 2,1kg</p>
	<p><b>Arm splint</b></p> <p><b>VA14X04</b> CLEAT with webbing 1-Chamber  <b>VA14X03</b> CLEAT with Webbing Multi-Chamber-Sack  <b>VA14CX04</b> CLICK with webbing 1-Chamber  <b>VA14CHX04</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VA14CX03</b> CLICK with Webbing Multi-Chamber-Sack  <b>VA14CHX03</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VA14CPX03</b> CLICK with Webbing DOT-Chamber  <b>VA14CPHX03</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 4 Width: top/hand: 50-30cm Length: 67cm Weight: 0,9kg</p>
	<p><b>Ankle splint</b></p> <p><b>VK15X04</b> CLEAT with webbing 1-Chamber  <b>VK15X03</b> CLEAT with Webbing Multi-Chamber-Sack  <b>VK15CX04</b> CLICK with webbing 1-Chamber  <b>VK15CHX04</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VK15CX03</b> CLICK with Webbing Multi-Chamber-Sack  <b>VK15CHX03</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VK15CPX03</b> CLICK with Webbing DOT-Chamber  <b>VK15CPHX03</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 4 Width: 47cm Length: 57cm Weight: 0,9kg</p>
	<p><b>Leg splint</b></p> <p><b>VL21X04</b> CLEAT with webbing 1-Chamber  <b>VL21X03</b> CLEAT with Webbing Multi-Chamber-Sack  <b>VL21CX04</b> CLICK with webbing 1-Chamber  <b>VL21CHX04</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VL21CX03</b> CLICK with Webbing Multi-Chamber-Sack  <b>VL21CHX03</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VL21CPX03</b> CLICK with Webbing DOT-Chamber  <b>VL21CPHX03</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 5 Width: top/hand: 70-40cm Length: 110cm Weight: 1,9kg</p>
	<p><b>Large leg splint</b></p> <p><b>VL18X04</b> CLEAT with webbing 1-Chamber  <b>VL18X03</b> CLEAT with Webbing Multi-Chamber-Sack  <b>VL18CX04</b> CLICK with webbing 1-Chamber  <b>VL18CHX04</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VL18CX03</b> CLICK with Webbing Multi-Chamber-Sack  <b>VL18CHX03</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VL18CPX03</b> CLICK with Webbing DOT-Chamber  <b>VL18CPHX03</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 5 Width: top/hand: 70-45cm Length: 135cm Weight: 2,1kg</p>







	<p><b>Forearm Splints</b></p> <p><b>VS100CKX02</b> Loop strap with CLICK 1-Chamber  <b>VS100CKX01</b> Loop strap with CLICK Multi-Chamber-Sack  <b>VS100CPKX01</b> Loop strap with CLICK DOT-Chamber</p> <p>Restraint straps: 2 Width: 35cm Length: 53cm Weight: 0,5kg</p>
	<p><b>Arm splints</b></p> <p><b>VS120CKX02</b> Loop strap with CLICK 1-Chamber  <b>VS120CKX01</b> Loop strap with CLICK Multi-Chamber-Sack  <b>VS120CPKX01</b> Loop strap with CLICK DOT-Chamber</p> <p>Restraint straps: 4 Width: 63-47cm Length: 74cm Weight: 1,1kg</p>
	<p><b>Leg Splints</b></p> <p><b>VS140CKX02</b> Loop strap with CLICK 1-Chamber  <b>VS140CKX01</b> Loop strap with CLICK Multi-Chamber-Sack  <b>VS140CPKX01</b> Loop strap with CLICK DOT-Chamber</p> <p>Restraint straps: 6 Width: 70-36cm Length: 110cm Weight: 2,1kg</p>
	<p><b>Forearm Splints</b></p> <p><b>VS100CX02</b> Webbing with CLICK 1-Chamber  <b>VS100CHX02</b> Wipe-Clean-Strap with CLICK 1-Chamber  <b>VS100CX01</b> Webbing with CLICK Multi-Chamber-Sack  <b>VS100CHX01</b> Wipe-Clean-Strap with CLICK Multi-Chamber-Sack  <b>VS100CPX01</b> Webbing with CLICK DOT-Chamber  <b>VS100CPHX01</b> Wipe-Clean-Strap with CLICK DOT-Chamber</p> <p>Restraint straps: 3 Width: 35cm Length: 53cm Weight: 0,5kg</p>
<p>Arm</p> 	<p><b>Arm splints</b></p> <p><b>VS120CX02</b> Webbing with CLICK 1-Chamber  <b>VS120CHX02</b> Wipe-Clean-Strap with CLICK 1-Chamber  <b>VS120CX01</b> Webbing with CLICK Multi-Chamber-Sack  <b>VS120CHX01</b> Wipe-Clean-Strap with CLICK Multi-Chamber-Sack  <b>VS120CPX01</b> Webbing with CLICK DOT-Chamber  <b>VS120CPHX01</b> Wipe-Clean-Strap with CLICK DOT-Chamber</p> <p>Restraint straps: 4 Width: 63-47cm Length: 74cm Weight: 1,1kg</p>
<p>Bein</p> 	<p><b>Leg Splints</b></p> <p><b>VS140CX02</b> Webbing with CLICK 1-Chamber  <b>VS140CHX02</b> Wipe-Clean-Strap with CLICK 1-Chamber  <b>VS140CX01</b> Webbing with CLICK Multi-Chamber-Sack  <b>VS140CHX01</b> Wipe-Clean-Strap with CLICK Multi-Chamber-Sack  <b>VS140CPX01</b> Webbing with CLICK DOT-Chamber  <b>VS140CPHX01</b> Wipe-Clean-Strap with CLICK DOT-Chamber</p> <p>Restraint straps: 6 Width: 70-36cm Length: 110cm Weight: 2,1kg</p>

	<p><b>Forearm Splints</b></p> <p><b>VS100CKX04</b> Loop strap with CLICK 1-Chamber  <b>VS100CKX03</b> Loop strap with CLICK Multi-Chamber-Sack  <b>VS100CPKX03</b> Loop strap with CLICK DOT-Chamber</p> <p>Restraint straps: 2 Width: 35cm Length: 53cm Weight: 0,5kg</p>
	<p><b>Arm Splints</b></p> <p><b>VS120CKX04</b> Loop strap with CLICK 1-Chamber  <b>VS120CKX03</b> Loop strap with CLICK Multi-Chamber-Sack  <b>VS120CPKX03</b> Loop strap with CLICK DOT-Chamber</p> <p>Restraint straps: 4 Width: 63-47cm Length: 74cm Weight: 1,1kg</p>
	<p><b>Leg Splints</b></p> <p><b>VS140CKX04</b> Loop strap with CLICK 1-Chamber  <b>VS140CKX03</b> Loop strap with CLICK Multi-Chamber-Sack  <b>VS140CPKX03</b> Loop strap with CLICK DOT-Chamber</p> <p>Restraint straps: 6 Width: 70-36cm Length: 110cm Weight: 2,1kg</p>
	<p><b>Forearm Splints</b></p> <p><b>VS100CX04</b> Webbing with CLICK 1-Chamber  <b>VS100CHX04</b> Wipe-Clean-Strap with CLICK 1-Chamber  <b>VS100CX03</b> Webbing with CLICK Multi-Chamber-Sack  <b>VS100CHX03</b> Wipe-Clean-Strap with CLICK Multi-Chamber-Sack  <b>VS100CPX03</b> Webbing with CLICK DOT-Chamber  <b>VS100CPHX03</b> Wipe-Clean-Strap with CLICK DOT-Chamber</p> <p>Restraint straps: 3 Width: 35cm Length: 53cm Weight: 0,5kg</p>
	<p><b>Arm Splints</b></p> <p><b>VS120CX04</b> Webbing with CLICK 1-Chamber  <b>VS120CHX04</b> Wipe-Clean-Strap with CLICK 1-Chamber  <b>VS120CX03</b> Webbing with CLICK Multi-Chamber-Sack  <b>VS120CHX03</b> Wipe-Clean-Strap with CLICK Multi-Chamber-Sack  <b>VS120CPX03</b> Webbing with CLICK DOT-Chamber  <b>VS120CPHX03</b> Wipe-Clean-Strap with CLICK DOT-Chamber</p> <p>Restraint straps: 4 Width: 63-47cm Length: 74cm Weight: 1,1kg</p>
	<p><b>Leg Splints</b></p> <p><b>VS140CX04</b> Webbing with CLICK 1-Chamber  <b>VS140CHX04</b> Wipe-Clean-Strap with CLICK 1-Chamber  <b>VS140CX03</b> Webbing with CLICK Multi-Chamber-Sack  <b>VS140CHX03</b> Wipe-Clean-Strap with CLICK Multi-Chamber-Sack  <b>VS140CPX03</b> Webbing with CLICK DOT-Chamber  <b>VS140CPHX03</b> Wipe-Clean-Strap with CLICK DOT-Chamber</p> <p>Restraint straps: 6 Width: 70-36cm Length: 110cm Weight: 2,1kg</p>




### Splints with Hook & Loop

	<p><b>Arm Splint</b></p> <p><b>VS101X01</b> Hook &amp; Loop straps with Multi-Chamber-Sack</p> <p>Restraint straps: 2            Width: Top/Hand: 40-30cm            Length: 73cm            Weight: 0,7kg</p>
	<p><b>Leg Splint</b></p> <p><b>VS103X01</b> Hook &amp; Loop straps with Multi-Chamber-Sack</p> <p>Restraint straps: 4            Width: Top/Foot: 60-18cm            Length: 96cm            Weight: 1,9kg</p>
	<p><b>Leg Splint</b></p> <p><b>VS105X01</b> Hook &amp; Loop straps with Multi-Chamber-Sack</p> <p>Restraint straps: 5            Width: Top/Foot: 73-18cm            Length: 133cm            Weight: 2,2kg</p>
	<p><b>Arm Splint</b></p> <p><b>VS101X02</b> Hook &amp; Loop straps with Multi-Chamber-Sack</p> <p>Restraint straps: 2            Width: Top/Hand: 40-30cm            Length: 73cm            Weight: 0,7kg</p>
	<p><b>Leg Splint</b></p> <p><b>VS103X02</b> Hook &amp; Loop straps with Multi-Chamber-Sack</p> <p>Restraint straps: 4            Width: Top/Foot: 60-18cm            Length: 96cm            Weight: 1,9kg</p>
	<p><b>Leg Splint</b></p> <p><b>VS105X02</b> Hook &amp; Loop straps with Multi-Chamber-Sack</p> <p>Restraint straps: 5            Width: Top/Foot: 73-18cm            Length: 133cm            Weight: 2,2kg</p>




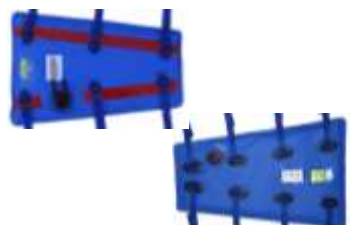


	<p><b>Forearm splint</b></p> <p><b>VS66X03</b> Cleats with Webbing 1-Chamber  <b>VS66CX03</b> CLICK with Webbing 1-Chamber  <b>VS66CHX03</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS67X03</b> Cleats with Webbing Multi-Chamber-Sack  <b>VS67CX03</b> CLICK with Webbing Multi-Chamber-Sack  <b>VS67CHX03</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VS67CPX03</b> CLICK with Webbing DOT-Chamber  <b>VS67CPHX03</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 3 Width: 30cm Length: 35cm Weight: 0.3ka</p>
	<p><b>Arm splint</b></p> <p><b>VS68X04</b> Cleats with Webbing 1-Chamber  <b>VS68CX04</b> CLICK with Webbing 1-Chamber  <b>VS68CHX04</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS69X03</b> Cleats with Webbing Multi-Chamber-Sack  <b>VS69CX03</b> CLICK with Webbing Multi-Chamber-Sack  <b>VS69CHX03</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VS69CPX03</b> CLICK with Webbing DOT-Chamber  <b>VS69CPHX03</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 4 Width: top/hand: 65-44cm Length: 70cm Weight: 1,5kg</p>
	<p><b>Leg splint</b></p> <p><b>VS70X03</b> Cleats with Webbing 1-Chamber  <b>VS70CX03</b> CLICK with Webbing 1-Chamber  <b>VS70CHX03</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS71X03</b> Cleats with Webbing Multi-Chamber-Sack  <b>VS71CX03</b> CLICK with Webbing Multi-Chamber-Sack  <b>VS71CHX03</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VS71CPX03</b> CLICK with Webbing DOT-Chamber  <b>VS71CPHX03</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 5 Width: top/foot: 70-55cm Length: 90cm Weight: 1,9kg</p>
	<p><b>Forearm splint</b></p> <p><b>VS66X04</b> Cleats with Webbing 1-Chamber  <b>VS66CX04</b> CLICK with Webbing 1-Chamber  <b>VS66CHX04</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS67X04</b> Cleats with Webbing Multi-Chamber-Sack  <b>VS67CX04</b> CLICK with Webbing Multi-Chamber-Sack  <b>VS67CHX04</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VS67CPX04</b> CLICK with Webbing DOT-Chamber  <b>VS67CPHX04</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 3 Width: 30cm Length: 35cm Weight: 0,3kg</p>
	<p><b>Arm splint</b></p> <p><b>VS68X05</b> Cleats with Webbing 1-Chamber  <b>VS68CX05</b> CLICK with Webbing 1-Chamber  <b>VS68CHX05</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS69X04</b> Cleats or Stripes with Webbing Multi-Chamber-Sack  <b>VS69CX04</b> CLICK with Webbing Multi-Chamber-Sack  <b>VS69CHX04</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VS69CPX04</b> CLICK with Webbing DOT-Chamber  <b>VS69CPHX04</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 4 Width: top/hand: 65-44cm Length: 70cm Weight: 1,5kg</p>
	<p><b>Leg splint</b></p> <p><b>VS70X04</b> Cleats with Webbing 1-Chamber  <b>VS70CX04</b> CLICK with Webbing 1-Chamber  <b>VS70CHX04</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS71X04</b> Cleats with Webbing Multi-Chamber-Sack  <b>VS71CX04</b> CLICK with Webbing Multi-Chamber-Sack  <b>VS71CHX04</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VS71CPX04</b> CLICK with Webbing DOT-Chamber  <b>VS71CPHX04</b> CLICK with Wipe-Clean Straps DOT-Chamber</p> <p>Restraint straps: 5 Width: top/foot: 70-55cm Length: 90cm Weight: 1,9kg</p>

## Vacuum-Mattresses

	<p><b>VM710X11</b> No patient restraint straps 1-Chamber  <b>VM710PX11</b> No patient restraint straps DOT-Chamber                  4 built in handles Width: 70cm Length: 100cm Weight: 3,5kg</p>
	<p><b>VM710X09</b> Cleats with Webbing 1-Chamber  <b>VM710X08</b> Cleats with Webbing Multi-Chamber-Sack  <b>VM710CX09</b> CLICK with Webbing 1-Chamber  <b>VM710CX08</b> CLICK with Webbing Multi-Chamber-Sack  <b>VM710CPX08</b> CLICK with Webbing DOT-Chamber  <b>VM710CHX09</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VM710CHX08</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VM710CPHX08</b> CLICK with Wipe-Clean Straps DOT-Chamber                  Carry handles: 4 Restraint straps: 3                  Width: 70cm Length: 100cm Weight: 3,5kg</p>
	<p><b>VM710X14</b> Cleats with Webbing 1-Chamber  <b>VM710X13</b> Cleats with Webbing Multi-Chamber-Sack  <b>VM710CX14</b> CLICK with Webbing 1-Chamber  <b>VM710CX13</b> CLICK with Webbing Multi-Chamber-Sack  <b>VM710CPX13</b> CLICK with Webbing DOT-Chamber  <b>VM710CHX14</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VM710CHX13</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VM710CPHX13</b> CLICK with Wipe-Clean Straps DOT-Chamber                  Carry handles: 4 Restraint straps: 3                  Width: 70cm Length: 100cm Weight: 4kg</p> <p><b>VM712X10</b> Cleats with Webbing 1-Chamber  <b>VM712X09</b> Cleats with Webbing Multi-Chamber-Sack  <b>VM712CX10</b> CLICK with Webbing 1-Chamber  <b>VM712CX09</b> CLICK with Webbing Multi-Chamber-Sack  <b>VM712CPX09</b> CLICK with Webbing DOT-Chamber  <b>VM712CHX10</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VM712CHX09</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VM712CPHX09</b> CLICK with Wipe-Clean Straps DOT-Chamber                  Carry handles: 4 Restraint straps: 3                  Width: 70cm Length: 120cm Weight: 4,5kg</p>

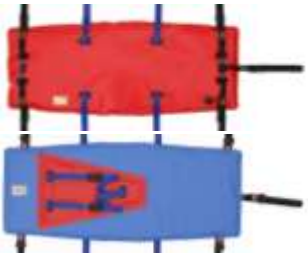

## Pediatric Vacuum Splints

	<p>Arm splint, angled  <b>VS56X20V</b> Hook &amp; Loop straps 1-Chamber                  Restraint straps: 3                  Width top/hand: 50-28cm Length: 42cm Weight: 0,6 kg</p>
	<p>Leg Splint  <b>VS58X20V</b> Hook &amp; Loop straps 1-Chamber                  Restraint straps: 3 Width top/foot: 38-25cm Length: 56cm Weight: 0,7kg</p>
	<p>Arm Splint, angled  <b>VS56X01</b> Cleats or Stripes with Webbing 1-Chamber  <b>VS56CX01</b> CLICK with Webbing 1-Chamber  <b>VS56CHX01</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS57CPX01</b> CLICK with Webbing DOT-Chamber  <b>VS57CPHX01</b> CLICK with Wipe-Clean Straps DOT-Chamber                  Restraint straps: 3                  Width top/hand: 50-28cm Length: 42cm Weight: 0,6 kg</p>
	<p>Leg Splint  <b>VS58X01</b> Cleats or Stripes with Webbing 1-Chamber  <b>VS58CX01</b> CLICK with Webbing 1-Chamber  <b>VS58CHX01</b> CLICK with Wipe-Clean Straps 1-Chamber  <b>VS59X01</b> Red Stripes with Webbing Multi-Chamber-Sack  <b>VS59CX01</b> CLICK with Webbing Multi-Chamber-Sack  <b>VS59CPX01</b> CLICK with Webbing DOT-Chamber  <b>VS59CPHX01</b> CLICK with Wipe-Clean Straps DOT-Chamber                  Restraint straps: 4 Width top/foot: 38-25cm Length: 56cm Weight: 0,7kg</p>

# RedVac Pediatric Systems

for use on the ambulance gurney

## Pediatric Vacuum-Mattresses




	<p> <b>VKF04X01</b> Cleats with Webbing Multi-Chamber-Sack  <b>VKF04CX01</b> CLICK with Webbing Multi-Chamber-Sack  <b>VKF04CHX01</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VKF04CPX01</b> CLICK with Webbing DOT-Chamber  <b>VKF04CPHX01</b> CLICK with Wipe-Clean Straps DOT-Chamber         </p> <p>           Restraint straps: 5    Securing straps: 5            Width Head/Middle/Foot: 53-61-53 cm    Length: 140cm    Weight: 3,5 kg         </p>
	<p> <b>VKF06X02</b> Cleats with Webbing Multi-Chamber-Sack  <b>VKF06CX02</b> CLICK with Webbing Multi-Chamber-Sack  <b>VKF06CHX02</b> CLICK with Wipe-Clean Straps Multi-Chamber-Sack  <b>VKF06CPX02</b> CLICK with Webbing DOT-Chamber  <b>VKF06CPHX02</b> CLICK with Wipe-Clean Straps DOT-Chamber         </p> <p>           Carry handles: 4    Restraint straps: 3            Width Top/Hand: 85-65 cm    Length: 115cm    Weight: 2,9 kg         </p>

## RedVac JUNIOR

Pediatric Restraint System for use with or without vacuum mattress



**JUNIOR** is a restraint system for safe transportation of children in ambulances. It is compatible with most ambulance gurneys. The system is suitable for new-borns and children up to around 7 years with body length 130 cm and weight of 35 kg.

<p><b>Junior</b></p> 	<p><b>KF01X07</b></p> <p>           Width head: 39 cm            Width foot: 20 cm            Length: 74 cm            Weight: 1,3 kg         </p> <p>           Restraint straps: 4            Securing straps: 6         </p>
<p>Vacuum mattresses for JUNIOR</p> 	<p> <b>VKF01PX02</b> Hook &amp; Loop straps DOT-Chamber  <b>VKF01CPX02</b> CLICK with Webbing DOT-Chamber  <b>VKF01CPHX02</b> CLICK with Wipe-Clean straps DOT-Chamber         </p> <p>           Width: 45 cm    Length: 80 cm    Weight: 0,75 kg         </p> <p>Restraint straps: 3</p>
	<p> <b>VKF02PX01</b> Hook &amp; Loop straps DOT-Chamber  <b>VKF02CPX01</b> CLICK with Webbing DOT-Chamber  <b>VKF02CPHX01</b> CLICK with Wipe-Clean straps DOT-Chamber         </p> <p>           Width Head/Middle/Foot: 45/62/62 cm    Length: 138 cm            Weight: 1,2 kg         </p> <p>Restraint straps: 5</p>

# Accessories for Vacuum Products

Can only be ordered together with the product

## Antiskid systems

**AS01**  
Three (3) stars with six cleats each



**AS02**  
Three (3) rubber bands



## Integrated Head Immobiliser

**HI01**



## Integrated Leg Separator

**LS01**



## Customised labels

Property indicator label



Label sizes:

8cm x 5cm

17cm x 7cm

A4

## Repair Kit

**9022**  
Repair Kit with  
self-adhesive patches



## Vacuum Pumps & Carry Bags

	<p><b>HP01</b> One-Hand pump, plastic Length: 21cm Width: 14cm Weight: 0,6kg</p>
	<p><b>HP03</b> Hand pump, aluminium Length: 30cm Width: 4cm Weight: 0,4kg</p>
	<p><b>HP04</b> Hand pump, plastic Length: 28cm Width: 11 / 4,5cm Weight: 0,170kg</p>
	<p><b>FP01</b> Foot pump, plastic Length: 43cm Width: 11cm Weight: 0,75kg</p>
	<p><b>HFP02</b> Hand- and Foot pump, plastic. Length: 45cm Width: 6cm Weight: 0,95kg</p>
	<p><b>FP04</b> Foot pump, aluminium Length: 54cm Width: 6cm Weight: 1,4kg</p>
<p>Carry bag for RedVac Junior</p> 	<p><b>TS04</b> Nylon fabric Colour: blue <b>TS04C B/R</b> Cordura fabric Colour: blue, red <b>TS04V B/R/G/Y</b> Wipe-off foil Colour: blue, red, green, yellow</p> <p>Length-Width-Height: 40-20-10cm Weight: 0,25kg</p>
<p>Carry bag for splint set</p> 	<p><b>TS06</b> Nylon fabric Colour: blue <b>TS06C B/R</b> Cordura fabric Colour: blue, red <b>TS06V B/R/G/Y</b> Wipe-off foil Colour: blue, red, green, yellow</p> <p>Length-Width-Height: 65-45-12cm Weight: 0,4kg</p>
<p>Carry bag for mattress with backpack straps</p>  	<p><b>TM01</b> Nylon fabric Colour: blue <b>TM01C B/R</b> Cordura fabric Colour: blue, red <b>TM01V B/R/G/Y</b> Wipe-off foil Colour: blue, red, green, yellow</p> <p>for mattress with width less than 90cm Length-Width-Height: 80-45-25cm Weight: 0,7kg</p> <p><b>TM02</b> Nylon fabric Colour: blue <b>TM02C B/R</b> Cordura fabric Colour: blue, red <b>TM02V B/R/G/Y</b> Wipe-off foil Colour: blue, red, green, yellow</p> <p>for mattress with width of more than 90cm Length-Width-Height: 100-45-30cm Weight: 1kg</p>

# TRS Thermo Rescue Sack

for use with a vacuum mattress  
and/or spine boards or rescue stretchers



- Heatinsulating, water repellent and rugged material provides good protection without additional blanket
- Designed for easy access to the patient
- 10 carry handles
- Under surface with antiskid stripes
- Zipper and Hook-and-Loop closing for fast handling and easy access to the patient
- Optional property indicator label
- CE marked



*Undersurface with antiskid stripes*




*Easy access to the patient*



*Also very tall patients can be accommodated*

## Specification

	<p><b>TRSX01</b> Material: Nylon Carry Handles: 10 Width Head/Middle/Foot: 45cm/70cm/40cm Length: 200cm Weight: 2,5kg</p>
	<p><b>TRSX02</b> Material: Nylon Carry Handles: 10 Width Head/Middle/Foot: 45cm/70cm/40cm Length: 220cm Weight: 2,9kg</p>
	<p><b>TRSX03</b> Material: CORDURA quilted Carry Handles: 10 Width Head/Middle/Foot: 45cm/70cm/40cm Length: 220cm Weight: 2,9kg</p>

# BRS Helicopter Rescue Sack

## for use with a vacuum mattress and/or spine boards or rescue stretchers



- Optimised for rescue with helicopter and use with vacuum mattress
- The vacuum mattress can be deflated without opening the rescue bag
- Antiskid under surface to prevent slipping on snow
- Hook and loop fastener and zipper for easy access to the patient and to avoid wind-sack phenomenon and rotation
- Stable head protection
- Lifting strap system is adjustable to achieve horizontal position of the patient
- Patient restraint straps with metal buckles
- Pocket for vacuum pump
- Plated-through holes for free drainage of water
- Anti-rotation rings
- CE labelled and certified for helicopter lifting

### BRS03

All straps are sewn onto the rescue bag and are stowed in easy to open pockets



### BRS04

The detachable lifting bridle is attached to the carry handles with karabiner



Easy access to the patient



Integrated bag for the vacuum pump



Undersurface with antiskid stripes



### Specification

	<p><b>BRS03</b>  <b>BRS04</b>          Carry Handles: 10          Patient restraint straps: 5          Width, Head-Middle-Foot: 40-70-45cm          Length: 220 cm          Weight: 4 kg</p> <p><b>TBRS03</b>          Carry bag          Width-Length-Height: 44-80-33 cm          Weight: 1 kg</p>
-------------------------------------------------------------------------------------	------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

# KOHLBRAT & BUNZ

***Emergency & Rescue Systems***

***Made in Austria***

Quality management in accordance with  
EN ISO 9001 - EN ISO 13485

certified by



*made by*

**KOHLBRAT & BUNZ GMBH**



**A-5550 Radstadt**

**Austria**

**T: +43 6452 7193 0**

**F: +43 6452 7193 51**

**office@kohlbrat-bunz.com**

**www.kohlbrat-bunz.com**

**www.redvac.com**

Pictures are symbolic and may differ from the original. Stated weights are approximate.  
Subject to technical changes.